

SPE100/200/400

Automatic Solid Phase Extraction

Shanghai Sineo SPE100/200/400 automatic solid phase extraction series adopts highly integrated pipeline design, with multi-function sample needle, relying on precise and flexible mechanical arm to improve sample preparation efficiency and bring you convenient and easy working experience, Let your work be faster.



Characteristics

The injection needle and the extraction column are sealed by a universal plug sealing structure and equipped with H-axis pressing block to ensure good sealing.

The suspension structure is adopted in the test tube bracket, which ensures the positioning accuracy of the manipulator to the greatest extent.

The instrument standard router can be connected by WiFi or network line, which fully realizes the whole process of automation and unattended.

Using modular design, only one communication line and pipeline can work between modules, and it is convenient to replace and maintain.

Add emergency stop button, when the instrument is abnormal, manual operation can be emergency stop, to maximize user safety. Handwheel adjusting origin greatly saves the time of adjusting origin.

Main features

Number of channels and control mode: The instrument host adopts modular design, supports the number of channels upgraded in the field, and supports 1/2/4/6/8 multiple channel combinations;

Continuous processing capacity:Maximum continuous processing capacity of 52-bit samples;

Robot arm:X/Y/Z three-dimensional manipulator design, X/Y axis linkage, improve operational efficiency; H axis is added to move the slider, giving consideration to flexible movement and accurate positioning;

Sample pump:Sampling accuracy, solvent corrosion resistance, high pressure resistance(maximum withstand voltage 1Mpa), and compatible with size and volume sampling;

Sample needle:316L stainless steel, solvent resistant; unique side-bottom opening to ensure full sample loading; highly integrated function for liquid channels, gas channels, and for extraction column seals;

Remote connection and control: The instrument is equipped with a router that can be connected via Wi-Fi or network cable, which fully realizes the whole process of automation and unattended operation;

Puncture sampling: Sample needle can be selected to realize puncture sampling and avoid sample volatilization loss;

Sample volume and receiving volume: The whole series supports 10 - 20mL sample tube, and can support 50mL and larger specifications within two channels; the receiving volume is standard 10mL/15mL, and the two channels can be expanded by 50mL;

Pressure sensing: Each channel is equipped with a standard pressure sensor, push column pressure to achieve closed-loop fuzzy control, and has the functions of overpressure alarm and termination of the experiment;

Velocity control: The top of each extraction column is equipped with an independent sealing cap. The high pressure sealing is realized by docking with the sample needle and moving Z-axis pressure block. The solution passes through the extraction column according to the set flow rate to ensure good reproducibility and parallelism;

Solvent management system: Supporting 8 reagents; Supporting the hint function of insufficient residual capacity; Supporting liquid level following function;

Infiltration and incubation: When the active solvent or eluting solvent is loaded into the extraction column, the liquid level can be kept above the packing for a certain time to ensure full activation or elution and to ensure the ideal experimental results;

Cleaning mode: Supporting various cleaning modes to ensure that the inner and outer walls of sample needles are clean without cross-contamination;

Waste discharge passage: Support environmental protection to recover waste liquor and separate treatment of waste liquor with different properties;

Specification of extraction column: Compatible with conventional 1/3/6mL extraction column, other specifications of extraction column can be customized according to requirements;

Multi-scheme extraction: Supporting the simultaneous operation of different specifications of extraction columns and different treatment projects without waiting, improving the efficiency of continuous treatment;

A/B column extraction: Support multi-mode tandem column extraction, automatic identification of A/B column extraction scheme;

Online drying/dewatering: Supporting on-line nitrogen purge mode;

Online Concentration: Supporting in-situ nitrogen purge and concentration, the concentration time can be edited;

Large volume injection: Supports the continuous and smooth loading of more than 500 mL samples without interruption, particle wear resistance, corrosion resistance, no need for any consumables and regular maintenance; supports the continuous switching of samples, the number of samples larger than 48 bits;

Cloud services: Supporting upload and download programs and parameters; Supporting online technical support and exchange; Supporting remote monitoring;

Booking starts: After editing the experimental scheme, the start-up time can be self-booked according to the time of the analyzer;

Count down: After the beginning of the experiment, the countdown time is used to facilitate the control of the end time of the experiment and to match the time of the analyzer in advance;


Intelligent self-check and alarm: After the parameter setting is completed, the software self-checks whether there are logic errors, promptly reminds the modification; Intelligent fault alarm system, according to the fault level, chooses to continue or stop;

Technical index

	SPE100	SPE200	SPE400	SPEX00 (Customization)
Channels	1	2	4	X
Channel Independent Control	Support			
Sample pump	Precision Continuous Sample Pump 0.1-120mL/min			
Pump flow rate	8Seed, 500 mL Solvent Bottle 0.1mL-20L			
Types of solvents	10 - 20 mL, customization			
Sample volume	10mL/15mL, other specifications customized Supporting			
Sample tube volume	X/Y/Z Three-Dimensional Manipulator 0.01mm			
Receiving volume	Docking Technology of Independent Seal Cover and Pan-plug Seal			
Puncture sampling	Positive pressure extraction			
Loading mode	80psi			
Accuracy of manipulator	Support, Overvoltage Alarm, Threshold Editable (5-145psi)			
Sealing technology	Nitrogen purge			
Extraction mode	Branch collection 3 channels Support, A/B Pillar			
column pressure	Nitrogen purge, solvent cleaning, blowing stirring cleaning and other modes are			
Pressure sensing	used to clean the inner and outer walls of sample needles.			
Column drying method	Support			
Waste collection	Solvent deficiency hints			
String string	680mm x 580mm x 670mm			
Cleaning mode	220/110V 50/60Hz 5-40Temperature			
	0-90% without condensation			
Liquid level following				
Solvent Residual Capacity				
Host Size (H*W*D)				
Voltage				
working temperature				
Working humidity				
power	60W	80W	120W	<200W

Your agent for Switzerland

Analytical instruments and services

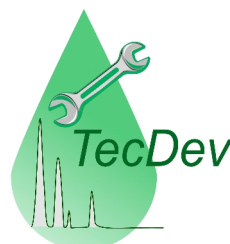

Le Villaret
Rte. de la Robellaz 6
CH - 1432 Belmont s/Yverdon
Switzerland

Tel: +41 (0) 24 435 22 02
Fax: +41 (0) 24 435 22 03
Mobile: +41 (0) 79 435 96 26
Email: tecdev@swissonline.ch
Visit us: www.tecdevexpress.com