

aquapuri 5 series

Water Purification System

The Most Needed WPS in Your Lab



CONTENTS

01

Types of
Water Grade

02

Advantages

03

Products

04

Selection Guide

01

Ttypes of Water Grade




Your agent for Switzerland:

TecDev
Le Villaret
Rte de la Robellaz 6
CH-1432 Belmont s/Yverdon
Switzerland

Phone: +41 (0) 24 435 22 02
Fax: +41 (0) 24 435 22 03
Mobile: +41 (0) 79 435 96 26
Email: tecdev@swissonline.ch
Visit us: www.tecdevexpress.ch

3 Types of Water

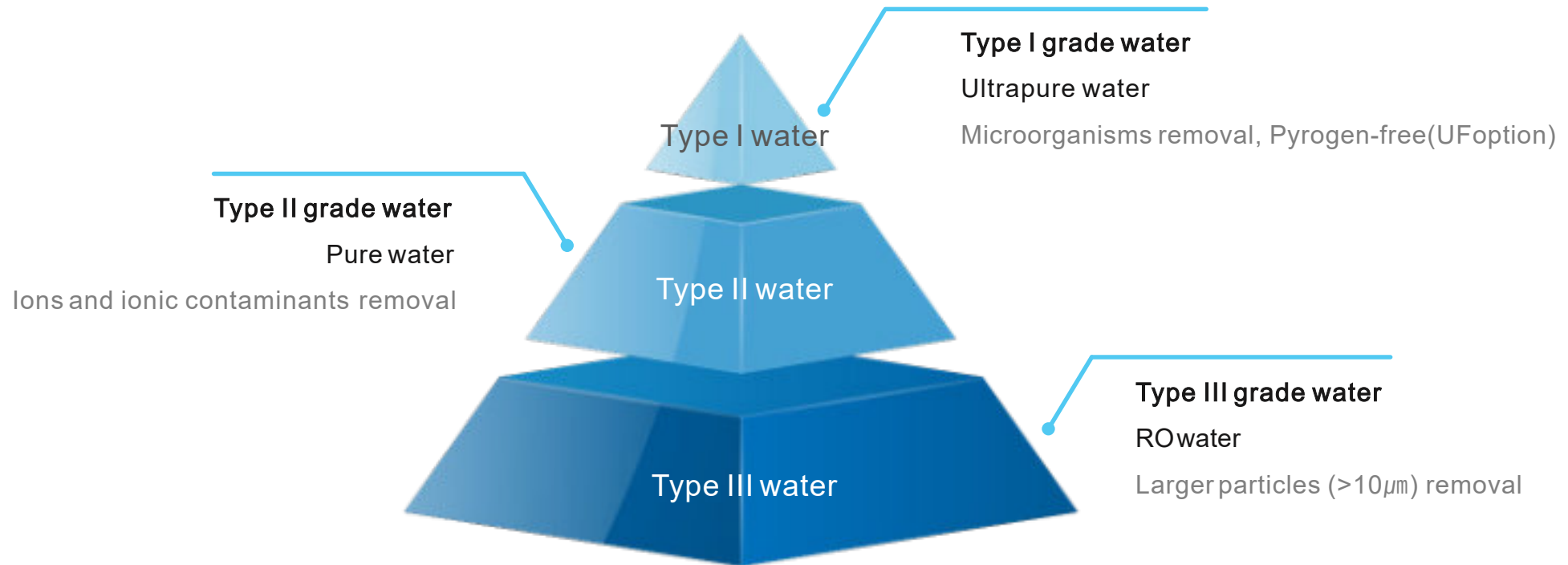

Key Factors of Water Purity

- Conductivity/Resistivity
- Organic Compound Levels
- Total Organic Carbon Levels
- Biological Contamination




Water Quality Parameters by ASTM

Parameter	Type I*	Type II**	Type III***
Conductivity ($\mu\text{S}/\text{cm}$) at 25°C, max	0.056	1.0	0.25
Resistivity ($\text{M}\Omega\text{-cm}$) at 25°C, max	18.0	1.0	4.0
pH at 25°C	—	—	—
TOC ($\mu\text{g}/\text{l}$), max	50	50	200
Sodium ($\mu\text{g}/\text{l}$), max	1	5	10
Silica ($\mu\text{g}/\text{l}$), max	3	3	500
Chloride ($\mu\text{g}/\text{l}$), max	1	5	10

*Requires use of 0.2 μm membrane filter; **Prepared by distillation; ***Requires the use of 0.45 μm membrane filter.

Flow Diagram of aquapuri System

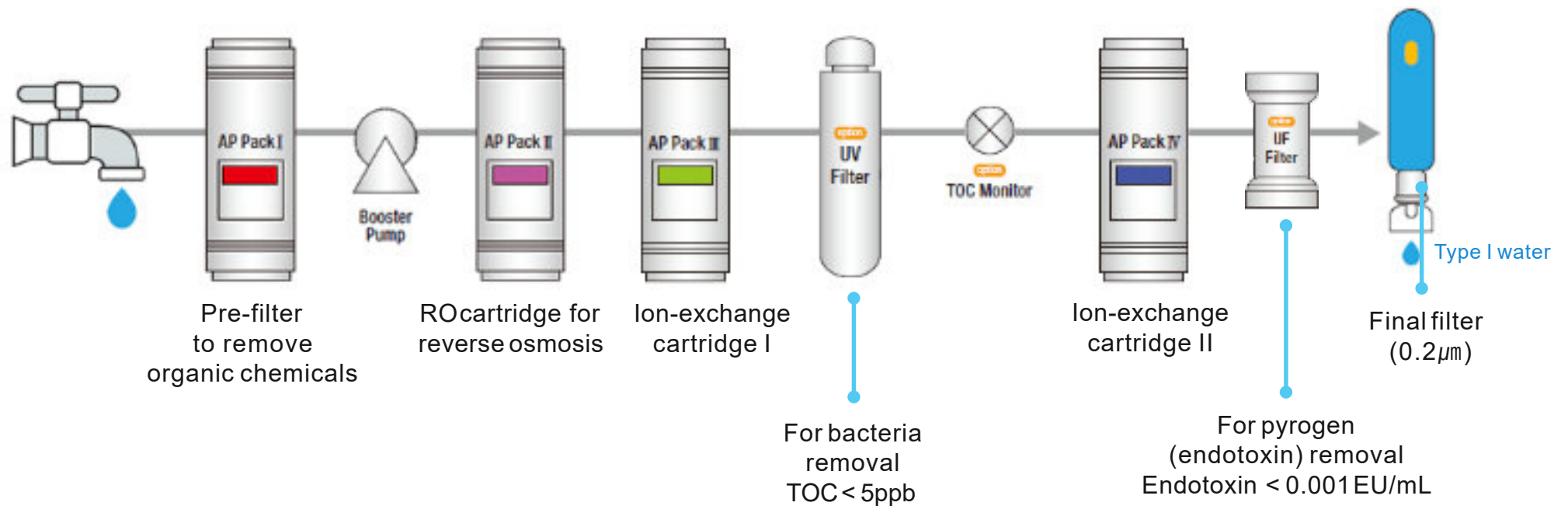

02

Advantages





Your agent for Switzerland:



TecDev
Le Villaret
Rte de la Robellaz 6
CH-1432 Belmont s/Yverdon
Switzerland

Phone: +41 (0) 24 435 22 02
Fax: +41 (0) 24 435 22 03
Mobile: +41 (0) 79 435 96 26
Email: tecdev@swissonline.ch
Visit us: www.tecdevexpress.ch

Advantages...No Tank, No Contamination

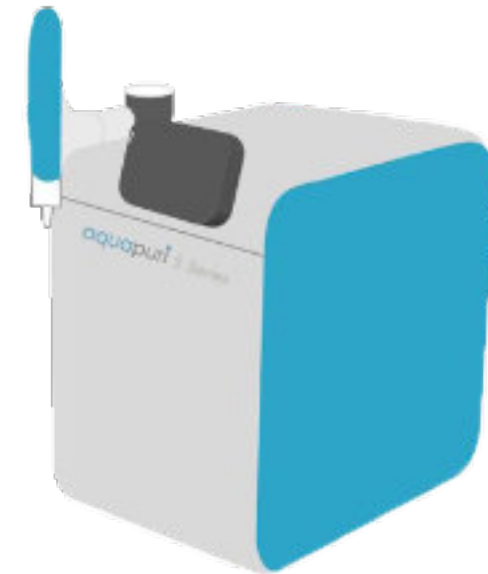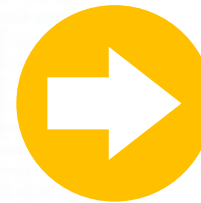


Advantages...Space Saving

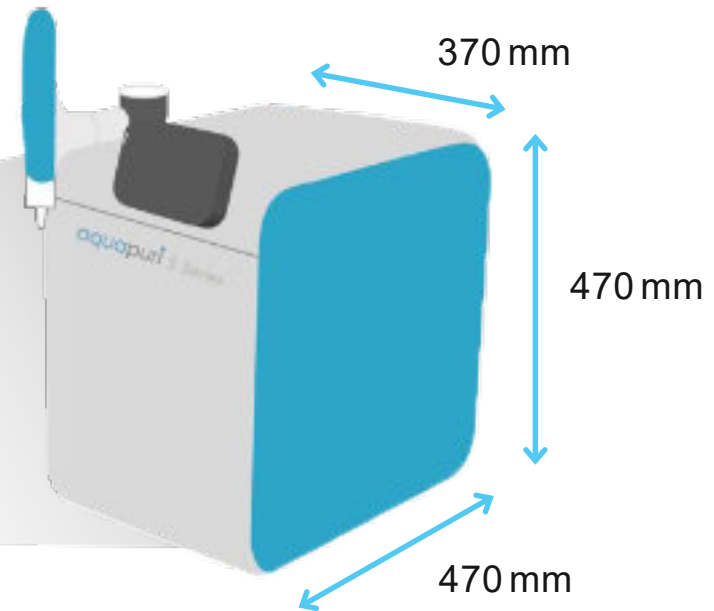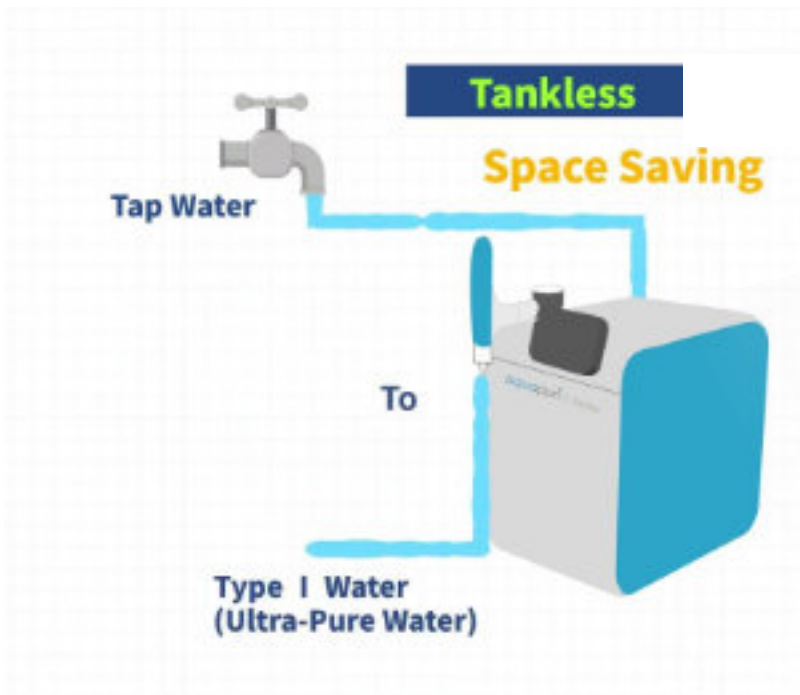




Advantages...Superior Productivity

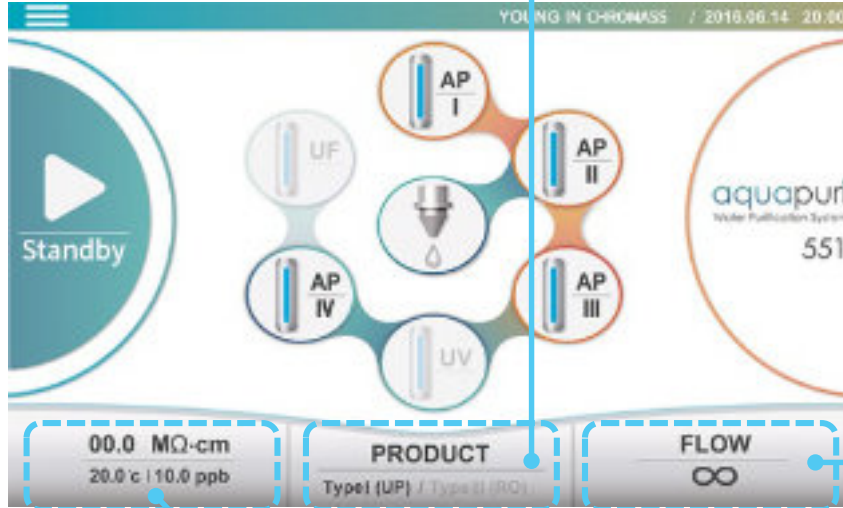


Water type selectable touch button

Flow rate setting
- Preset volume
- Manual volume

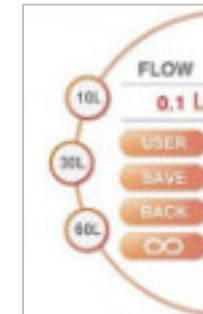

Water quality, temperature and TOC level monitoring (Switchable to monitor resistivity/conductivity)

Advantages...Intuitive Touch LCD

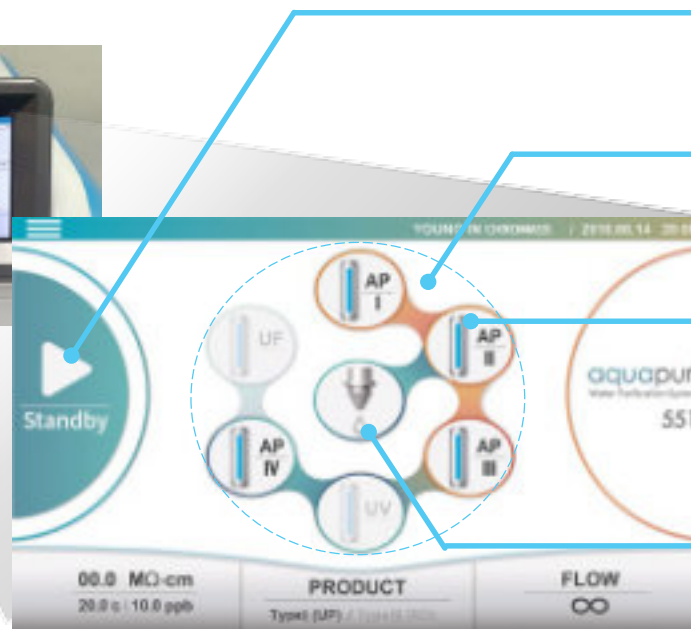

Intuitive touch-screen interface

Standby/Operate button

Replacement alert
/ Auto-flush for cleaning

Life cycle monitoring on filters



All filters are changeable at once

Advantages...User Convenience





Advantages...Tool-Free Maintenance

Easy door open



Magnetic door

Filter Replacement without Tools



03

Products



Your agent for Switzerland:

TecDev
Le Villaret
Rte de la Robellaz 6
CH-1432 Belmont s/Yverdon
Switzerland

Phone: +41 (0) 24 435 22 02
Fax: +41 (0) 24 435 22 03
Mobile: +41 (0) 79 435 96 26
Email: tecdev@swissonline.ch
Visit us: www.tecdevexpress.ch

Type II



363 series



521 series

Type I



551 series



532 series

Type III



541 series



511 series

aquapuri 5 series



aquapuri 363

Features

Type II water

- Type II grade water production
- Resistivity: $10 \text{ M}\Omega \cdot \text{cm}$ at 25°C
- Conductivity: $0.1 \mu\text{s}/\text{cm}$ at 25°C
- Delivery Flow Rate: 20 L/h



aquapuri 511/521

Features of aquapuri 511

Type III water

- Conductivity: $10 \mu\text{s}/\text{cm}$ at 25°C
- Delivery Flow Rate: 60 L/h

Features of aquapuri 521

Type II water

- Resistivity: $10 \text{ M}\Omega \cdot \text{cm}$ at 25°C
- Conductivity: $0.1 \mu\text{s}/\text{cm}$ at 25°C
- Delivery Flow Rate: 60 L/h



aquapuri 532

Features

Type I water

- Resistivity: $18.2 \text{ M}\Omega \cdot \text{cm}$ at 25°C
- Conductivity: $0.055 \mu\text{S}/\text{cm}$ at 25°C
- TOC: $< 5 \text{ ppb}$ (Option)
- Endotoxin (Pyrogen): $< 0.001 \text{ EU}/\text{mL}$ (Option)
- Delivery Flow Rate: $120 \text{ L}/\text{h}$
- Type II water required as feed water



aquapuri 541

Features

Integrated model (Type I, III)

- Tank-less system
- Resistivity: $18.2 \text{ M}\Omega \cdot \text{cm}$ at 25°C
- Conductivity: $0.055 \mu\text{s/cm}$ at 25°C
- TOC: $< 5 \text{ ppb}$ (Option)
- Endotoxin (Pyrogen): $< 0.001 \text{ EU/mL}$ (Option)
- Delivery Flow Rate: 60L/h



aquapuri 551

Features

Integrated model (Type I, II)

- Tank-less system
- Resistivity: $18.2 \text{ M}\Omega \cdot \text{cm}$ at 25°C
- Conductivity: $0.055 \mu\text{S}/\text{cm}$ at 25°C
- TOC: $< 5 \text{ ppb}$ (Option)
- Endotoxin (Pyrogen): $< 0.001 \text{ EU}/\text{mL}$ (Option)
- Delivery Flow Rate: $60\text{L}/\text{h}$

04

Selection Guide



Your agent for Switzerland:

TecDev
Le Villaret
Rte de la Robellaz 6
CH-1432 Belmont s/Yverdon
Switzerland

Phone: +41 (0) 24 435 22 02
Fax: +41 (0) 24 435 22 03
Mobile: +41 (0) 79 435 96 26
Email: tecdev@swissonline.ch
Visit us: www.tecdevexpress.ch

- Washing/Rinsing
- Autoclaves
- Steam Generators

- Buffer/Media prep.
- Electrophoresis
- Feed water for ultra-pure water

- Reagents
- General analysis (IC, AAS, HPLC, ICP-MS)
- Standard solvents

- Ultra-trace inorganic/organic analysis
- GC-MS and TOC analysis

- Mammalian cell culture
- DNA sequencing
- Monoclonal antibody production

- PCR
- 2-Dimensional electrophoresis
- Cell culture

Type I

532
541
551

532UV
541UV
551UV

532UF
541UF
551UF

532VF
541VF
551VF

Type II

363
521
551(Type II)

Type III

511
541(Type III)

